

Multiple socket outlet e-medical™ Q400AN EU




Premium aluminium multiple socket outlet with overload protection

Rock solid, square, excellent: The MED4-Q400AN EU multiple socket is a real space-saving marvel and impresses with quality and safety in a wide range of medical usage. For example, for use directly at patient's bedsides, in operating theatres or on mobile cart solutions for endoscopy or mobile rounds. 8 fuses protect valuable equipment from disruptions in the electrical circuit. In addition, the MED4-Q400AN EU has an extra overload fuse. This ensures, that the main switch does not trip, even if the total load of 3500 W is exceeded. The MED4-Q400AN EU multiple socket provides one PE connector pin per consumer and can be used as a PE distribution.

Benefits:

- **Square 4-times** multiple socket outlet
- Schuko sockets
- **2 x fuse** per socket insert
- Supply cables with a cross-section of 2.5 mm²
- TÜV tested for use with anaesthesia equipment


Features

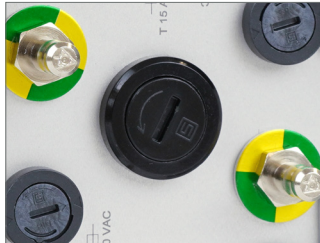

1x additional fuse against overload



2 fuses per socket insert. Replaceable during operation.

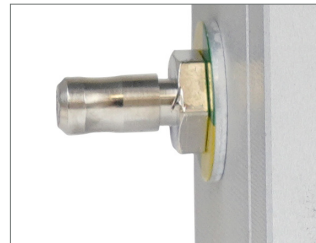

Additional potential equalization (PE) on the housing



11 mm thick power cable with 2.5 mm² single lines

Technical data

e-medical™ MED4-Q400AN EU

Number of Plugs	4
Connector type	Schuko
Fuse per socket insert	2 x 4A microfuses inserts each 5x 20 mm (DIN EN 60127-2/VDE 0820-2)
Additional fuse	1x 15A microfuse inserts 5x 20 mm (DIN 41.571-3)
ZPA (PE) connector	according to DIN 42801
Load per socket	max. 850 W
Total load multiple outlet strip	max. 3600W
Current max.	15 A
Input voltage	230 V
Dimensions (LWH)	190 x 150 x 83 mm
Case	Aluminium
Standards	IEC 60601-1:2005 + A1:2012 IEC 60601-1:2005 + A1:2012/AMD1:2012 ISO 80601-2-13:2011/AMD1:2005
Warranty	5 years
Scope of delivery	1x multiple socket outlet, 2x MED4 Double guard bracket, 1x socket wrench
3 m cable length	Item No: 2012098
5 m cable length	Item No: 2012142